



Accesories for Operating Tables

FENIX OT-01 VIVAX OT-02



ACCESORIES FOR OPERATING TABLES

This catalogue presents the wide range of accessories available for INFIMED operating table.

They were created basing on long experience in designing of medical equipment, hours spent on consultation with medical staff and functional and ergonomy testing.

The accessories haven been sorted according to practical use and surgical disciplines.


Table of contents:



General surgery positions



SUPINE POSITION

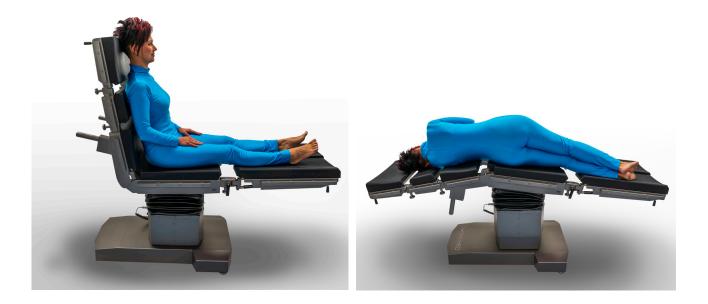

BEACH CHAIR POSITION

LATERAL POSITION AND FLEX





EB-01 One position clamp

TECHNICAL INFORMATION:

For fixation of the various accesories with rods of 16 mm diameter to the side rail

Made of stainless steel, with a matt finish



EB-02 Multi position clamp

TECHNICAL INFORMATION:

For fixation of the various accesories with rods of 16 mm diameter to the side rail

Allows radial adjustment of the accesories

Made of stainless steel, with a matt finish





EB-04 Angle arm support

TECHNICAL INFORMATION:

Height adjustment: 300mm

Distance from table top adjustment: 200mm

Tilting adjustment by ball joint

Made of stainless steel with matt finish and polyurethane foam, with two fixing velcro straps

To be fixed with EB-02 clamp

EB-05 Arm support

TECHNICAL INFORMATION:

Height adjustment: 100mm

Made of stainless steel with matt finish and polyurethane foam, with two fixing velcro straps

To be fixed with EB-01 clamp





TECHNICAL INFORMATION: Height adjustment: 300mm Made of stainless steel with matt finish and polyurethane foam To be fixed with EB-01 clamp



EB-07 Side-shoulder support (side support position)

TECHNICAL INFORMATION: Height adjustment: 230mm Distance from table top adjustment: 150mm Made of stainless steel with matt finish and polyurethane foam To be fixed with EB-01 clamp

EB-07 Side-shoulder support (shoulder support position)

TECHNICAL INFORMATION: Height adjustment: 230mm Distance from table top adjustment: 150mm Made of stainless steel with matt finish and polyurethane foam To be fixed with EB-01 clamp



EB-08 Side support

TECHNICAL INFORMATION: Height adjustment: 150mm Distance from table top adjustment: 150mm Made of stainless steel with matt finish and polyurethane foam To be fixed with EB-01 clamp









EB-09 Anesthetic frame

TECHNICAL INFORMATION: Height adjustment: 600mm Made of stainless steel with matt finish To be fixed with EB-01 clamp

EB-10 Anesthetic frame with adjustable width

TECHNICAL INFORMATION:

Height adjustment: 500mm

Both side witdth adjustment 200mm Made of stainless steel with matt finish

To be fixed with EB-01 clamp



EB-11 Infusion holder

TECHNICAL INFORMATION:

Height adjustment: 200mm

Made of stainless steel with matt finish

To be fixed with EB-01 clamp



EB-12 Hand (shank holder)

TECHNICAL INFORMATION:

Height adjustment: 300mm

Made of stainless steel with matt finish and polyurethane foam

To be fixed with EB-01 clamp



EB-13 Support for meniscus operation

TECHNICAL INFORMATION: Height adjustment: 600mm Made of stainless steel with matt finish and polyurethane foam To be fixed with EB-02 clamp



EB-14 Arm support, with horizontal movement

TECHNICAL INFORMATION:

Height adjustment: 300mm

Distance from tabletop adjustment: 200mm

Tilting adjustment by ball joint

Made of stainless steel with matt finish, HPL, foam mattress covered with antistatic artificial leather, fixing velcro straps To be used with EB-01 clamp

ES-01 Leg strap ES-02 Thigh strap

TECHNICAL INFORMATION:

Adjustment by velcro: 400mm

Strap width: 100mm

Made of stainless steel with matt finish and polyamide strap



ES-03 Abdominal strap

TECHNICAL INFORMATION:

Adjustment by velcro: 400mm

Strap width: 100mm

Made of stainless steel with matt finish and polyamide strap







ES-04 Hand strap

TECHNICAL INFORMATION: Adjustment by velcro: 400mm Strap width: 100mm

Made of stainless steel with matt finish and polyamide strap

ES-07 Holder for additional equipment, mobile

TECHNICAL INFORMATION:

Non-marking wheels with brakes

Hooks and rails for assembly of accessories on rack

Basket in basis

Made of stainless steel with matt finish

ES-08 Holder for arm operations

TECHNICAL INFORMATION:

Fixed to side rails of operating table

Vertical adjustment of angle against tabletop Distance from tabletop adjustment: 100mm

Height adjustable feet

X-Ray permeable tabletop, dimensions 800x400mm

Made of stainless steel with matt finish, foam mattress covered with antistatic artificial leather

ES-09 Holder for tubes

TECHNICAL INFORMATION:

Height adjustment: 200mm

Distance from tabletop adjustment: 200mm

Made of stainless steel with matt finish and HPL

To be fixed with EB-01 clamp







ES-10 Tray for additional accessories

TECHNICAL INFORMATION: Height adjustment: 400mm Tray dimensions: 450x250mm Made of stainless steel with matt finish To be used with EB-01 clamp



ES-11 X-ray tray

TECHNICAL INFORMATION: Height adjustment: 400mm To be used with X-Ray cassette of width381mm Installation from head or side Made of stainless steel with matt finish



ES-12 Chest support

TECHNICAL INFORMATION: Height adjustment: 150mm Distance from tabletop adjustment: 150mm Made of stainless steel with matt finish and polyurethane foam To be fixed with EB-01 clamp



ES-13 Attachement for knee operation

TECHNICAL INFORMATION: Height adjustment: 200 mm Knee holder adjustment Made of stainless steel with matt finish and polyurethane foam To be fixed with EB-01 clamp





ES-14 Pubic support

TECHNICAL INFORMATION: Height adjustment: 150 mm Distance from tabletop adjustment: 150 mm Made of stainless steel with matt finish and polyurethane foam To be fixed with EB-01 clamp



ES-15 Leg support, not-divided

TECHNICAL INFORMATION:

To be fixed to siderails

Made of stainless steel with matt finish and foam mattress covered with antistatic artificial leather



ES-17 Side rail extension

TECHNICAL INFORMATION:

Extension for side rail for mounting side rail accessories

Made of stainless steel with matt finish



ES-18 Tabletop extension

TECHNICAL INFORMATION:

To be fixed to siderails

Made of stainless steel with matt finish, HPL and polyure than foam



ES-19 Strap for anesthetic frame

TECHNICAL INFORMATION: To be fixed to anesthetic frame EB-09 or EB-10 Made of polymide strap



ES-20 Side extension for bariatric tabletop

TECHNICAL INFORMATION:

To be fixed to siderails

Made of stainless steel with matt finish, HPL and polyurethane foam with artificial letter



ES-21 Leg support, divided

TECHNICAL INFORMATION:

To be fixed to siderails

Made of stainless steel with matt finish and polyurethane foam



Gynaecology Urology Proctology Positions



GYNEACOLOGICAL SURGERY



PROCTOLOGICAL (KNEEL) SURGERY

Gynaecology Urology Proctology Accesories



EB-03 Knee support

TECHNICAL INFORMATION: Height adjustment: 100mm Made of stainless steel with matt finish and polyurethane foam To be fixed with EB-02 clamp



ES-05 Gyneacological bowl

TECHNICAL INFORMATION: Fixed instead of the leg supports Adjustable horizontal position and rotation Bowl capacity: 4 litres

Made of stainless steel with matt finish



ES-06 Urological bowl

TECHNICAL INFORMATION: Fixed instead of the leg supports Adjustable horizontal position and rotation Bowl capacity:4 litres, with rubber outflow Made of stainless steel with matt finish



ES-16 Proctological attachement

TECHNICAL INFORMATION:

Height adjustment: 200mm

Kneeling chair, without clamp, fixed to side rails

Abdominal support to be used with 2 pcs of EB-02 clamps

Made of stainless steel with matt finish, HPL and polyure than foam



Neurosurgery Laryngology Ophthalmology Positions



NEUROSURGICAL SURGERY - SUPINE POSITION (WITH DORO NEUROSURGICAL ATTACHEMENT)



NEUROSURGICAL SURGERY - SITTING POSITION (WITH DORO NEUROSURGICAL ATTACHEMENT)

Neurosurgery Laryngology Ophthalmology Accesories



EA-01 Adapter for specialistic head rest

TECHNICAL INFORMATION:

Neccesary to fix to the operating table to use DORO neurosurgical attachement or any of specialistic head rest

Made of stainless steel with matt finish



EN-01 Specialistic head rest with adjustable track

TECHNICAL INFORMATION:

Fixed to the adapter EA-01

Height adjustment 150mm, angle adjustment to the table top

Adjustment of the distance from the tabletop: 250mm

Made of stainless steel with matt finish and polyurethane foam





EN-02 Specialistic head rest, horseshoe

TECHNICAL INFORMATION: Fixed to the adapter EA-01 Height adjustment 150mm, angle adjustment to the table top Adjustment of the distance from the tabletop: 250mm Made of stainless steel with matt finish and polyurethane foam

EN-05 Arm support

TECHNICAL INFORMATION:

Height adjustment: 300mm

Distance from tabletop adjustment: 200mm

Tilting adjustment by ball joint

Made of stainless steel with matt finish, HPL, foam mattress covered with antistatic artificial leather, fixing velcro straps To be used with EB-01 clamp





DORO neurosurgical attachement for laying position

TECHNICAL INFORMATION:

To be fixed to adapter EA-01

Including:

- Adjustable Base Unit 3001-00 •
- Swivel Adaptor 3002-00
- QR3 Skull Clamp 1001.001
- Reusable Skull Pins Pediatric 3004-00 (3 pcs) •
- Reusable Skull Pins Adult 3005-00 (3 pcs)

DORO neurosurgical attachement for sitting position


TECHNICAL INFORMATION:

To be fixed to side rails using EB-02 multiposition clamps

Including:

- Adjustable Base Unit 3001-00
- Swivel Adaptor 3002-00 QR3 Skull Clamp 1001.001
- Reusable Skull Pins Pediatric 3004-00 (3 pcs)
- Reusable Skull Pins Adult 3005-00 (3 pcs)
- Crossbar Adaptor

Orthopaedics Traumatology Positions

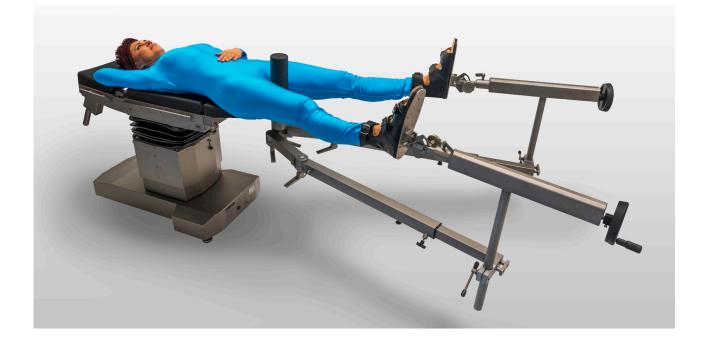

ORTHOPAEDIC SURGERY - SUPINE POSITION (ORTHOPAEDIC ATTACHEMENT EO-01)



ORTHOPAEDIC SURGERY - SUPINE POSITION (ORTHOPAEDIC ATTACHEMENT EO-01 WITH SUPPORTS EO-10 AND BAR SUPPORT EO-02)



Orthopaedics Traumatology Positions



ORTHOPAEDIC SURGERY - LATERAL DECUBITUS POSITION (ORTHOPAEDIC ATTACHEMENT EO-01 WITH SUPPORTS EO-10 AND HIP SUPPORT EO-09)

Ortopaedic attachemnt EO-01 (with EO-10 supports)

TECHNICAL INFORMATION: Includes 2 extension devices with foot plates and straps Adidtional accessories (E0-02...EO-11) are to be ordered seperately Fixed to leg rests sockets Adjustment of extensions devices height: 350mm Distance between support roll and extensions devices: 300-1300mm Adjustment of rotation arms position: 0-180 Adjustment of length of otrhopaedic attachments arms: 600mm Support roll X-Ray permeable Made of stainless steel with matt finish, foam mattress covered with artificial leather, polyurethane

Orthopaedics Traumatology Accesories



EO-02 Bar support

TECHNICAL INFORMATION:

To be fixed to the arms of the orthopeadic attachement EO-01 with EO-07 clamp

Height adjustment 400mm

Made of stainless steel with matt finish and leather strap



EO-03 Leg support

TECHNICAL INFORMATION:

To be fixed to the arms of the orthopeadic attachement EO-01 with EO-07 clamp

Height adjustment 400mm

Made of stainless steel with matt finish and leather strap / polyure-thane foam



EO-04 Support for plastering

TECHNICAL INFORMATION:

To be fixed to the socket of the orthopeadic attachement EO-01

Made of stainless steel with matt finish and leather strap / polyure-thane foam



EO-05 Heel holder

TECHNICAL INFORMATION:

To be fixed to orthopeadic attachement EO-01 instead of plate supports

Made of stainless steel with matt finish







EO-06 foot support

TECHNICAL INFORMATION:

To be fixed to the orthopeadic attachement EO-01 instead of plate supports

Made of stainless steel with matt finish and leather straps

EO-07 Clamp for orthopaedic attachment supports

TECHNICAL INFORMATION:

Used to fix the following accesories: EO-02, EO-03, EO-08, EO-11, to the orthopeadic attachement EO-01

Made of stainless steel with matt finish





EO-08 Knee operation attachment for orthopedic attachment

TECHNICAL INFORMATION: Without extension device To be fixed to the arms of the orthopeadic attachement EO-01 with EO-07 clamp Height adjustment: 200mm Adjustment of distance between support roll and extension device: 800mm Angle adjustment against tabletop Made of stainless steel with matt finish and polyurethane foam

EO-09 Hip support for lateral decubitus position

TECHNICAL INFORMATION:

To be fixed to the socket of the orthopaedic attachment EO-01

Adjustment left/right

Adjustment of support roll rotation

Made of stainless steel with matt finish and polyurethane foam

Orthopaedics Traumatology Accesories



EO-10 Support for orthopedic attachment arms

TECHNICAL INFORMATION:

To be fixed to the arms of the orthopeadic attachement EO-01

Height adjustment: 400 mm

Made of stainless steel with matt finish

EO-11 Knee support for side positions

TECHNICAL INFORMATION:

To be fixed to the arms of the orthopeadic attachement EO-01 with EO-07 clamp

Made of stainless steel with matt finish and polyurethane foam





TECHNICAL INFORMATION:

Used to transport and assembly the orthopeadic attachement EO-01 to the operating table

Non-marking castors with individual brakes

Basket for accesories in the base

Made of stainless steel with matt finish

EO-13 Mobile stand for additional equipment of orthopedic attachement

TECHNICAL INFORMATION:

Non-marking castors with individual brakes

Hooks and rails for storage of orthopaedic accessories

Basket in the basi

Made of stainless steel with matt finish and HPL





Orthopaedics Traumatology Accesories



EO-14 Leather shoes for orthopedic attachment

TECHNICAL INFORMATION:

To be fixed to the orthopeadic attachement EO-01 instead of plate supports

Made of stainless steel with matt finish and leather



EO-20 Attachment for hand operating

TECHNICAL INFORMATION:

To be fixed with 2 pcs of clamp EB-01

Height adjustment: 200mm Adjustment of extension device length: 350mm

Adjsutment of tilt and rotation of hand fixing

Made of stainless steel with matt finish and polyurethane foam

EO-30 Attachement for shoulder operations head rest EN-02

TECHNICAL INFORMATION:

To be fixed to legrests sockets

Includes EN-02 head rest Removable right and left shoulder sections

Made of stainless steel with matt finish, polyurethane foam and HPL



Angiography table



VIEW FROM THE SIDE, WITH MAX LONGITUDINAL SLIDE



VIEW FROM THE BACKSIDE, WITH MAX LATERAL SLIDE





Angiography Accesories

ER-01.0

TECHNICAL INFORMATION: To be fixed on carbo fibre tabletop Allows to fix heavier additional accesories on the side rails Made of stainless steel with matt finish.



TECHNICAL INFORMATION: To be fixed on carbo fibre tabletop Allows to fix lighter additional accesories on the side rails Made of stainless steel with matt finish.



ER-02.0

TECHNICAL INFORMATION:

To be fixed on carbo fibre tabletop

Allows to fix heavier additional accesories on the side rails

X-ray permeable

Made of carbon fibre and stainless steel with matt finish.



ER-02.1

TECHNICAL INFORMATION:

To be fixed on carbo fibre tabletop

llows to fix lighter additional accesories on the side rails

X-ray permeable

Made of carbon fibre and stainless steel with matt finish.

Positioning pads

EM-01 Mattress for spine surgeries

TECHNICAL INFORMATION:

For positioning the patient, especially during surgical procedures on the spine

Made of polyurethane foam with artificial letter



EM-02 Mattress for sitting position

TECHNICAL INFORMATION:

For positioning the patient, especially during procedures requiring sitting position


Made of polyurethane foam with artificial letter



Gel patient body positioners

TECHNICAL INFORMATION:

Thick layered gel patient body positioners for patient on operation table: arm, chest, heel, head positioners and table pads

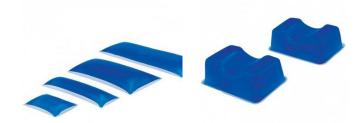





European Union

European Regional Development Fund



For special requirement of the customer it is possible to produce the product with changed technical parameters which do not diminish its safety.

The details given in this catalogue are correct at the time of going to press. Infimed Sp. z o.o., however, reserves the right to improve shown equipment.

adres INFIMED Sp. z o.o. ul. Kabaty 1, 34-300 Żywiec, Polska tel/fax +48 33 861 40 96 email office@infimed.pl www www.infimed.pl
 NIP
 553-251-29-67

 REGON
 243274947

 KRS
 0000463933

Sąd Rejonowy w Bielsku-Białej, VIII Wydział Gospodarczy Krajowego Rejestru Sądowego. Wysokość Kapitału Zakładowego: 500 000 PLN bank Bank Spółdzielczy w Węgierskiej Górce SWIFT POLUPLPR

78 8131 0005 0016 2492 2000 0010 PLN PL02 8131 0005 0016 2492 2000 0020 EUR PL23 8131 0005 0016 2492 2000 0030 USD







INFIMED 2014/C-Tables-EN-v1.2